




Rabbit Anti-Goat IgG(H+L), Human SP ads

Cat. No.	Format	Size
6164-01	Purified (UNLB)	1.0 mg
6164-02	Fluorescein (FITC)	1.0 mg
6164-04	Alkaline Phosphatase (AP)	1.0 mL
6164-05	Horseradish Peroxidase (HRP)	1.0 mL
6164-08	Biotin (BIOT)	1.0 mg

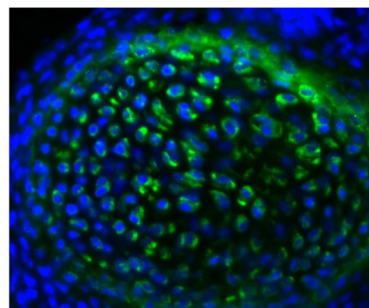

Frozen newborn mouse rib section was stained with Goat Anti-Type II Collagen-UNLB (SB Cat. No. 1320-01) followed by Rabbit Anti-Goat IgG(H+L), Human SP ads-FITC (SB Cat. No. 6164-02) and DAPI.

Description

Specificity	Reacts with the heavy and light chains of goat IgG
Source	Pooled antisera from rabbits hyperimmunized with goat IgG
Cross Adsorption	Human serum proteins (SP); may react with immunoglobulins from other species and the light chains of other goat immunoglobulins
Purification	Affinity chromatography on goat IgG covalently linked to agarose

Applications

Quality tested applications include –

ELISA¹
FLISA

Other referenced applications include –

FC²⁻⁶
IHC-FS⁷
IHC-PS⁸⁻¹⁰
ICC¹¹
WB¹²⁻¹⁴

Working Dilutions

ELISA	AP conjugate	1:2,000 – 1:4,000
	HRP conjugate	1:4,000 – 1:8,000
	BIOT conjugate	1:5,000 – 1:20,000

FLISA	FITC conjugate	1:200 – 1:400
--------------	----------------	---------------

Other Applications Since applications vary, you should determine the optimum working dilution for the product that is appropriate for your specific need.

For Research Use Only. Not for Diagnostic or Therapeutic Use.

Corporate Offices: 160 Oxmoor Blvd • Birmingham, AL 35209 • USA Mailing Address: P.O. Box 26221 • Birmingham, AL 35260 • USA

Tel: 205.945.1774 • U.S. and Canada: 800.722.2255 • Fax: 205.945.8768

Email: info@southernbiotech.com • Website: www.southernbiotech.com

Handling and Storage

- The purified (UNLB) antibody is supplied as 1.0 mg purified immunoglobulin in 1.0 mL of borate buffered saline, pH 8.2. *No preservatives or amine-containing buffer salts added.* Store at 2-8°C.
- The fluorescein (FITC) conjugate is supplied as 1.0 mg in 1.0 mL of PBS/NaN₃. Store at 2-8°C.
- The alkaline phosphatase (AP) conjugate is supplied as 1.0 mL in a stock solution of 50 mM Tris/1 mM MgCl₂/50% glycerol, pH 8.0, containing NaN₃ as preservative. Store at 2-8°C or long-term at -20°C.
- The horseradish peroxidase (HRP) conjugate is supplied as 1.0 mL in a stock solution of 50% glycerol/50% PBS, pH 7.4. No preservative added. Store at 2-8°C or long-term at -20°C.
- The biotin (BIOT) conjugate is supplied as 1.0 mg in 2.0 mL of PBS/NaN₃. Store at 2-8°C.
- Protect fluorochrome-conjugated forms from light. Reagents are stable for the period shown on the label if stored as directed.

Warning

Some reagents contain sodium azide. Please refer to product specific SDS.

References

1. Boxx GM, Kozel TR, Nishiya CT, Zhang MX. Influence of mannan and glucan on complement activation and C3 binding by *Candida albicans*. *Infect Immun.* 2010;78:1250-9. (ELISA)
2. Fu S, Gong F, Xie C, Zhu W, Wang W, Shen H, et al. S100A12 on circulating endothelial cells surface in children with Kawasaki disease. *Pediatr Res.* 2010;68:165-8. (FC)
3. Rezaee F, Rellick SL, Piedimonte G, Akers SM, O'Leary HA, Martin K, et al. Neurotrophins regulate bone marrow stromal cell IL-6 expression through the MAPK pathway. *PLoS One.* 2010;5(3):e9690. (FC)
4. Sekiya I, Ojima M, Suzuki S, Yamaga M, Horie M, Koga H, et al. Human mesenchymal stem cells in synovial fluid increase in the knee with degenerated cartilage and osteoarthritis. *J Orthop Res.* 2012;30:943-9. (FC)
5. Gong F, Zhang Y, Xie C, Zhu W, Wang W, Fu S, et al. Expression of receptor for advanced glycation end products (RAGE) on the surface of circulating endothelial cells is upregulated in Kawasaki disease. *Pediatr Res.* 2012;71:720-4. (FC)
6. Hegde M, Corder A, Chow KK, Mukherjee M, Ashoori A, Kew Y, et al. Combinational targeting offsets antigen escape and enhances effector functions of adoptively transferred T cells in glioblastoma. *Mol Ther.* 2013;21:2087-101. (FC)
7. Kreft KL, Verbraak E, Wierenga-Wolf AF, van Meurs M, Oostra BA, Laman JD, et al. The IL-7R α pathway is quantitatively and functionally altered in CD8 T cells in multiple sclerosis. *J Immunol.* 2012;188:1874-83. (IHC-FS)
8. Crijns AP, de Graeff P, Geerts D, ten Hoor KA, Hollema H, van der Sluis T, et al. MEIS and PBX homeobox proteins in ovarian cancer. *Eur J Cancer.* 2007;43:2495-505. (IHC-PS)
9. Ando K, Brion J, Stygelbout V, Suain V, Authélet M, Dedecker R, et al. Clathrin adaptor CALM/PICALM is associated with neurofibrillary tangles and is cleaved in Alzheimer's brains. *Acta Neuropathol.* 2013;125:861-78. (IHC-PS)
10. Ahmed L, Nalwoga H, Arnes JB, Wabinga H, Micklem DR, Aklsen LA. Increased tumor cell expression of Axl is a marker of aggressive features in breast cancer among African women. *APMIS.* 2015 May 25. doi: 10.1111/apm.12403. [Epub ahead of print]. (IHC-PS)
11. Hosper NA, van den Berg PP, de Rond S, Popa ER, Wilmer MJ, Masereeuw R, et al. Epithelial-to-mesenchymal transition in fibrosis: collagen type I expression is highly upregulated after EMT, but does not contribute to collagen deposition. *Exp Cell Res.* 2013;319:3000-9. (ICC)
12. Rayes J, Roumenina LT, Dimitrov JD, Repessé Y, Ing M, Christophe O, et al. The interaction between factor H and VWF increases factor H cofactor activity and regulates VWF prothrombotic status. *Blood.* 2014;123:121-5. (WB)
13. Wu X, Sagave J, Rutkovskiy A, Haugen F, Baysa A, Nygård S, et al. Expression of bone morphogenetic protein 4 and its receptors in the remodeling heart. *Life Sci.* 2014;97:145-54. (WB)
14. Matsebatiela TM, Anderson AL, Gallicchio VS, Elford H, Rice CD. 3,4-Dihydroxy-benzohydroxamic acid (Didox) suppresses pro-inflammatory profiles and oxidative stress in TLR4-activated RAW264.7 murine macrophages. *Chem Biol Interact.* 2015;233:95-105. (WB)