
# SouthernBiotech 1





# Goat F(ab')<sub>2</sub> Anti-Human Ig

Cat. No.	Format	Size
2012-01	Purified (UNLB)	0.5 mg
2012-02	Fluorescein (FITC)	0.5 mg
2012-04	Alkaline Phosphatase (AP)	1.0 mL
2012-05	Horseradish Peroxidase (HRP)	1.0 mL
2012-08	Biotin (BIOT)	0.5 mg
2012-09	R-phycoerythrin (PE)	0.25 mg
2012-30	Alexa Fluor® 488 (AF488)	0.5 mg
2012-32	Alexa Fluor® 555 (AF555)	0.5 mg

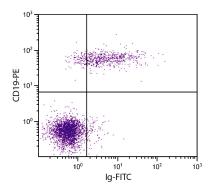

Human peripheral blood lymphocytes were stained with Goat F(ab')<sub>2</sub> Anti-Human Ig-FITC (SB Cat. 2012-02) and Mouse Anti-Human CD19-PE (SB Cat. No. 9340-09).

## **Description**

Specificity Reacts with the heavy and light chains of human IgG, IgM, and IgA

Source Pepsin digest of Goat Anti-Human Ig (SB Cat. No. 2010)

Cross Adsorption None; may react with immunoglobulins from other species

# **Applications**

Quality tested applications include -

ELISA 1,2,5 FLISA FC 7,8

Other referenced applications include -

ELISPOT <sup>3,4</sup> IHC-FS <sup>6</sup> Stim <sup>9,10</sup>

### **Working Dilutions**

ELISA	AP conjugate	1:1,000 – 1:4,000
	HRP conjugate	1:2,000 - 1:8,000
	BIOT conjugate	1:5,000 - 1:20,000
FLISA	FITC, AF488, and AF555 conjugates	1:100 – 1:400
	PE conjugate	$\leq$ 1 $\mu$ g/mL
Flow Cytometry	FITC, BIOT, and AF488 conjugates	$\leq$ 1 $\mu$ g/10 $^6$ cells
	PE conjugate	$\leq 0.1 \ \mu g/10^6 \ cells$
	For flow cytometry, the suggested use of these reagents is in a final volume of 100 $\mu\text{L}$	
Other Applications	Since applications vary, you should determine the optimum working dilution for the product that is appropriate for your specific need.	

For Research Use Only. Not for Diagnostic or Therapeutic Use.

#### **Handling and Storage**

- The purified (UNLB) antibody is supplied as 0.5 mg purified immunoglobulin in 1.0 mL of borate buffered saline, pH 8.2. No
  preservatives or amine-containing buffer salts added. Store at 2-8°C.
- The fluorescein (FITC), Alexa Fluor® 488 (AF488), and Alexa Fluor® 555 (AF555) conjugates are supplied as 0.5 mg in 1.0 mL of PBS/NaN<sub>3</sub>. Store at 2-8°C.
- The alkaline phosphatase (AP) conjugate is supplied as 1.0 mL in a stock solution of 50 mM Tris/1 mM MgCl₂/50% glycerol, pH 8.0, containing NaN₃ as preservative. Store at 2-8°C or long-term at -20°C.
- The horseradish peroxidase (HRP) conjugate is supplied as 1.0 mL in a stock solution of 50% glycerol/50% PBS, pH 7.4. No preservative added. Store at 2-8°C or long-term at -20°C.
- The biotin (BIOT) conjugate is supplied as 0.5 mg in 1.0 mL of PBS/NaN<sub>3</sub>. Store at 2-8°C.
- The R-phycoerythrin (PE) conjugate is supplied as 0.25 mg in 1.0 mL of PBS/NaN₃ and a stabilizing agent. Store at 2-8°C. Do not freeze!
- Protect fluorochrome-conjugated forms from light. Reagents are stable for the period shown on the label if stored as directed.

#### Warning

Some reagents contain sodium azide. Please refer to product specific SDS.

#### References

- Kozlowski PA, Chen D, Eldridge JH, Jackson S. Contrasting IgA and IgG neutralization capacities and responses to HIV type 1 gp120 V3 loop in HIV-infected individuals. AIDS Res Hum Retroviruses. 1994;10:813-22. (ELISA)
- Sydor JR, Normant E, Pien CS, Porter JR, Ge J, Grenier L, et al. Development of 17-allylamino-17-demethoxygeldanamycin hydroquinone hydrochloride (IPI-504), an anti-cancer agent directed against Hsp90. Proc Natl Acad Sci USA. 2006;103:17408-13. (ELISA)
- 3. Sira MM, Yoshida T, Takeuchi M, Kashiwayama Y, Futatani T, Kanegane H, et al. A novel immunoregulatory protein in human colostrum, syntenin-1, for promoting the development of IgA-producing cells from cord blood B cells. Int Immunol. 2009;21:1013-23. (ELISPOT)
- Rehnberg M, Amu S, Tarkowski A, Bokarewa MI, Brisslert M. Short- and long-term effects of anti-CD20 treatment on B cell ontogeny in bone marrow
  of patients with rheumatoid arthritis. Arthritis Res Ther. 2009;11(4):R123. (ELISPOT)
- 5. White CA, Pone EJ, Lam T, Tat C, Hayama KL, Li G, et al. Histone deacetylase inhibitors upregulate B cell microRNAs that silence AID and Blimp-1 expression for epigenetic modulation of antibody and autoantibody responses. J Immunol. 2014;193:5933-50. (ELISA)
- Xu W, He B, Chiu A, Chadburn A, Shan M, Buldys M, et al. Epithelial cells trigger frontline immunoglobulin class switching through a pathway regulated by the inhibitor SLPI. Nat Immunol. 2007;8:294-303. (IHC-FS)
- 7. Heile JM, Fong Y, Rosa D, Berger K, Saletti G, Campagnoli S, et al. Evaluation of hepatitis C virus glycoprotein E2 for vaccine design: an endoplasmic reticulum-retained recombinant protein is superior to secreted recombinant protein and DNA-based vaccine candidates. J Virol. 2000;74:6885-92. (FC)
- 8. Bao G, Li Y, Ma Q, He X, Xing J, Yang X, et al. Isolating human antibody against human hepatocellular carcinoma by guided-selection. Cancer Biol Ther. 2005;4:1374-80. (FC)
- 9. Liubchenko GA, Appleberry HC, Striebich CC, Franklin KE, Derber LA, Holers VM, et al. Rheumatoid arthritis is associated with signaling alterations in naturally occurring autoreactive B-lymphocytes. J Autoimmun. 2013;40:111-21. (Stim)
- 10. Lyubchenko T, Zerbe GO. B cell receptor signaling-based index as a biomarker for the loss of peripheral immune tolerance in autoreactive B cells in rheumatoid arthritis. PLoS One. 2014;9(7):e102128. (Stim)

Alexa Fluor® 488, 647, and 555 are provided under an Intellectual property license from Life Technologies Corporation. The purchase of this product conveys to the buyer the non-transferable right to use the purchased amount of the product and components of the product in research conducted by the buyer (whether the buyer is an academic or for-profit entity). The buyer cannot sell or otherwise transfer (a) this product (b) its components or (c) materials made using this product or its components to a third party or otherwise use this product or its components or materials made using this product or its components for Commercial Purposes. Commercial Purposes means any activity by a party for consideration and may include, but is not limited to: (1) use of the product or its components in manufacturing; (2) use of the product or its components to provide a service, information, or data; (3) use of the product or its components for therapeutic, diagnostic or prophylactic purposes; or (4) resale of the product or its components, whether or not such product or its components are resold for use in research. For information on purchasing a license to this product for any other use, contact Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, CA 92008 USA or outlicensing@lifetech.com.

Email: info@southernbiotech.com • Website: www.southernbiotech.com