













Horizontal Electrophoresis - CE Series (BL-Mini-CE)

www.biolinkk.com


























Horizontal Electrophoresis System (BL-Mini-CE)

The Biolinkk Horizontal Electrophoresis System (Model BL-Mini-CE) is a compact and efficient unit designed for small to medium DNA and RNA gel electrophoresis, featuring UV transparent gel trays, adjustable gel levels, and easy maintenance.



- UV Transparent Gel Tray.
- Adjustable Gel Level
- Replaceable Electrode Cassettes
- Lock-in Design
- Leak-free casting
- Impact, High Temperature & corrosion resistance

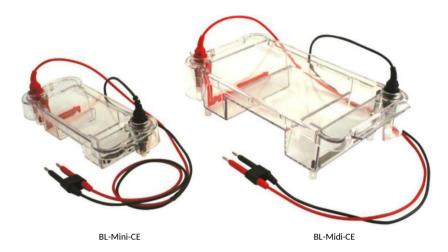

The Biolinkk Horizontal Electrophoresis System, Model BL-Mini-CE, is ideal for laboratories requiring precise and reliable DNA and RNA gel electrophoresis. This system includes a UV transparent gel tray with a fluorescent ruler for easy visualization, and the gel tray has level adjustability to suit various experiment needs. With a sample output capacity of 8-30 samples and a buffer volume of 270 ml, the system ensures efficient separation of nucleic acids. The clear PC molding allows easy observation of samples during the run. It also features easy-to-replace electrode cassettes and a lock-in design for running precast Agarose gels, making it user-friendly and highly functional.



- UV Transparent Gel Tray: Includes a fluorescent ruler for easy sample measurement.
- Adjustable Gel Level: Provides flexibility to suit different experimental requirements.
- Sample Output: Capable of processing 8-30 samples per run.
- Buffer Volume: 270 ml to maintain stable pH during electrophoresis.
- Clear PC Molding: Ensures clear visibility of samples during electrophoresis.
- Replaceable Electrode Cassettes: Facilitates easy maintenance and quick replacement.
- Lock-in Design: Designed for easy use with precast Agarose gels.







Specifications	BL-Mini-CE
Gel Size (cm) (W X L)	7 X 7, 7 X 10
Gel Tray	Yes
Gel Caster	Level Adjustable
Sample Output	8 - 30
Buffer Volume (ml)	270
Cell Size (cm)	26X12X6.5
Distance to Electrode (cm)	14.5



Add: Third Floor, Plot No. 25, Community Centre, Block B-1, Janakpuri, Delhi-110058

To order contact us:

Email: customerservice@biolinkk.com Tel: +91-11-20003-5753 / 25553658 www.bolinkk.com





