




Horizontal Electrophoresis System (BL-SPHT)

www.biolinkk.com


 Innovation to support science....

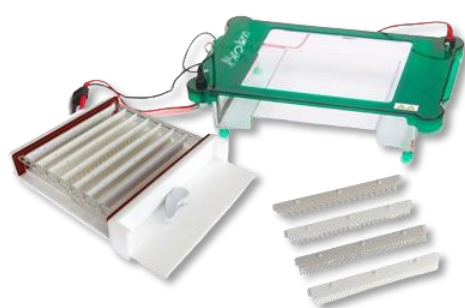


Horizontal Electrophoresis System (BL-SPHT)

The Biolinkk Horizontal Electrophoresis System (Model BL-SPHT) is a compact and efficient unit designed for laboratory experiments. It ensures stable pH, easy maintenance, and safety, making it ideal for DNA and RNA gel electrophoresis.



Highlights

- ◆ Compact 245×250mm gel size for efficient nucleic acid separation.
- ◆ Stable pH and cooling effects with a 2000 ml buffer volume.
- ◆ Transparent lid with auto-switch-off for safety.
- ◆ Easy maintenance with unfoldable electrodes and color-coded tray.

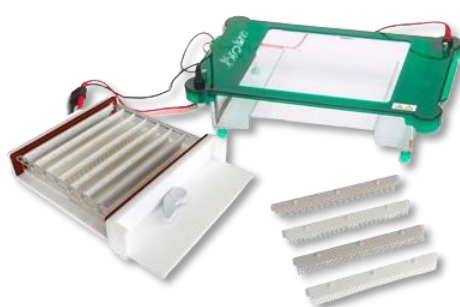

The Biolinkk Horizontal Electrophoresis System, Model BL-SPHT, offers precision and reliability in nucleic acid separation. With a gel size of 245×250mm and a buffer volume of 2000 ml, it provides a consistent cooling effect and stable pH throughout the experiment. The system features a transparent top lid to prevent liquid volatilization and ensures safe operation with an auto-switch-off mechanism when the lid is removed. The unfolded electrode design allows for easy repairs and part replacements, while the colored tray background facilitates the convenient addition of samples. This system is ideal for routine DNA and RNA gel electrophoresis applications in molecular biology labs.



Features

- ◆ Gel Size: 245×250mm.
- ◆ Buffer Volume: 2000 ml for optimal cooling and pH stability.
- ◆ Transparent top lid to prevent liquid evaporation and electricity leakage.
- ◆ Auto-switch-off safety feature when lid is removed.
- ◆ Unfolded electrodes for easy repair and part replacement.
- ◆ Sample background color design for ease in adding samples.
- ◆ Adding sample quicker and more suitable for the multi-channel pipettes.





Specifications

Cat. No.	BL-SPHT
Gel Size (W×L)	245×250mm
Sample volume	19~364
Buffer Volume	2000ml
Dimension (L×W×H) mm	487×300×93
Net weight	3.8 (kg)



Add: Third Floor, Plot No. 25, Community Centre, Block B-1, Janakpuri, Delhi-110058

To order contact us:

Email: customerservice@biolinkk.com
Tel: +91-11-20003-5753 / 25553658
www.biolinkk.com

